

# Q6 edge HD

Get the Edge.



**Q6 Edge® HD** shown with  
TRU-Balance® 3 HD Power  
Positioning and Q-Logic 2  
Drive Controls

With features like standard 4-pole motors, Mid-Wheel 6® Drive Design and ATX Suspension, the Q6 Edge® HD is engineered to meet the performance needs of the most active user. The Q6 Edge HD accepts a complete range of seating and electronics options to deliver exceptional rehab capability at an excellent value.



**QUANTUM®**

The Rehab Power Chair Solutions Company

[www.quantumrehab.com](http://www.quantumrehab.com)


# Specifications

<b>DRIVE WHEELS</b>	35.56 cm (14") solid (standard) or pneumatic (optional)
<b>CASTOR WHEELS</b>	<b>Front:</b> 15.24 cm (6") solid <b>Rear:</b> 15.24 cm (6") solid
<b>SUSPENSION</b>	ATX Suspension (Active-Trac® with extra stability) with enhanced rear suspension
<b>MAXIMUM SPEED<sup>1</sup></b>	Up to 7.6 km/h (4.7 mph)
<b>BRAKES</b>	Intelligent braking (electronic, regenerative disc brakes)
<b>GROUND CLEARANCE<sup>2</sup></b>	7.62 cm (3") at base frame
<b>TURNING RADIUS</b>	50.2 cm (19.75")
<b>SIZE<sup>2</sup></b>	Length: 90.4 cm (35.6") (without foot riggings) Width: 66.68 cm (26.25") (without foot riggings)
<b>SEATING</b>	Synergy® Seat Adjustable Width: 51 cm - 71 cm (20"-28") Adjustable Depth: 51 cm - 71 cm (20"-28") TRU-Balance® 2 HD Power Positioning Adjustable Width: 51 cm - 71 cm (20"-28") Adjustable Depth: 51 cm - 71 cm (20"-28") TRU-Balance® 3 HD Power Positioning Adjustable Width: 51 cm - 71 cm (20"-28") Adjustable Depth: 51 cm - 63.5 cm (20"-25") High-back Seat Width/Depth: 51 cm - 71 cm (20"-28") Medium-back Seat Width/Depth: 51 cm - 71 cm (20"-28")
<b>DRIVETRAIN</b>	Mid-Wheel 6®, 4-pole
<b>BATTERIES</b>	70 AH (Group 24)
<b>RANGE PER CHARGE<sup>1,3</sup></b>	Up to 25 km (15.5 miles)
<b>BATTERY CHARGER</b>	8A, off-board
<b>AVAILABLE ELECTRONICS</b>	100A Q-Logic Drive Control System
<b>WEIGHT CAPACITY</b>	204 kg (450 lbs.)
<b>BASE WEIGHT</b>	65 kg (143.5 lbs.) (without batteries)
<b>BATTERY WEIGHT<sup>4</sup></b>	23.6 kg (52 lbs.) (each)
<b>WARRANTY</b>	Lifetime limited warranty on frame; 2-year warranty on electronics; 18-month warranty on drive motor

# Shroud Colours



Candy Apple Red    Viper Blue    Onyx Black    Sunburst Orange    True Blue

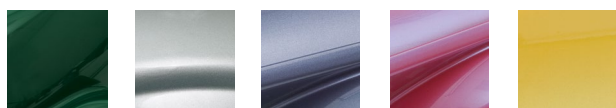

Emerald Green    Silver    Sapphire Blue    Pearl Pink    Orion Yellow

# Features

- Mid-Wheel 6® allows six wheels on the ground for maximum stability
- Compatible with TRU-Balance® 3 HD Power Positioning
- Low impact OMNI-Castors (nylon, spherical-shaped castors) on front and rear prevent wheel hang-ups
- ATX Suspension (Active-Trac® with extra stability) incorporates front OMNI-Castors and semi-independent rear castor beam for enhanced performance over more varied terrain

# Options

- TRU-Balance® 3 (450 lbs.) HD Tilt, HD Recline, and HD Tilt & Recline
- Headrest
- Occupied transit for TRU-Balance® 2 HD Tilt <sup>5</sup>
- Occupied transit for Synergy® Seating <sup>5</sup>

### FOOTNOTES:

- 1) Range and speed vary with user weight, terrain type, battery charge, battery condition and tyre pressure.
- 2) Due to manufacturing tolerances and continual product improvement, this specification can be subject to variance of + or - 3%. Dimensions listed are for power base only. Overall measurements will vary based on seating and accessory selections.
- 3) Tested in accordance with ANSI/RESNA, WC Vol. 2, Section 4, standards. Results derived from theoretical calculations based on battery specifications and drive system performance. Testing conducted at maximum weight capacities. This specification can be subject to a variance of +10%, -5%.
- 4) Battery weight may vary based on manufacturer.
- 5) Limited to 24" W x 24" D.

The information contained herein is correct at the time of publication; due to Pride's commitment to constant improvement and development, we reserve the right to alter specifications without prior notice.

Theoretical range is calculated under test conditions in accordance with European standards. Actual range in normal use will depend on many factors, including the condition of the vehicle and its batteries, the weight of the driver, correct tyre pressures, ambient temperature, and the gradient and surface of the road or pavement.



— LIVE YOUR BEST! —

AU\_Q6 Edge HD\_SS\_7-09-14

20-24 Apollo Drive  
Hallam, Victoria 3803  
www.pridemobility.com.au

Tel: 03 8770 9600  
Fax: 03 9703 2720



The Rehab Power Chair Solutions Company



Quality Endorsed Company  
8001 QEC20026  
SAI Global