


# Dedicated to Revolutionising the Power Chair Experience







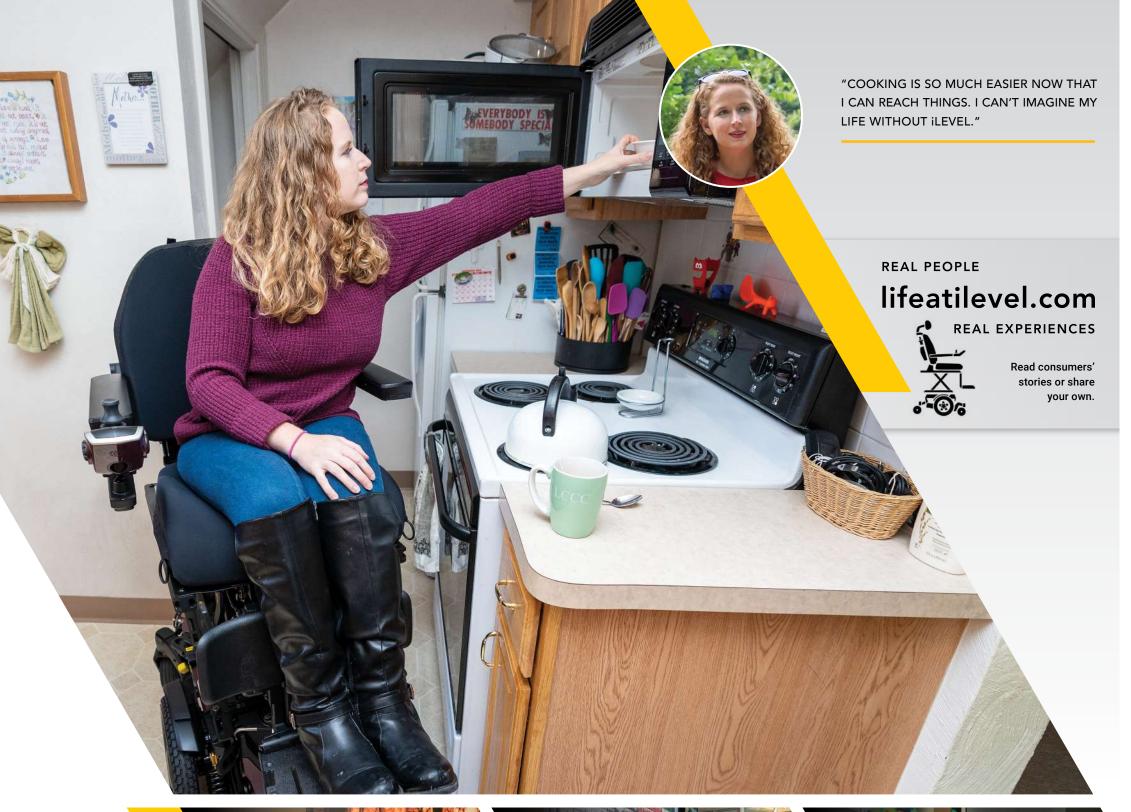

# Safer transfers Improved reach interaction at eye level #1 FOR REHAB POWER

# iLEVEL®

# **POWER SEAT ELEVATION**

# Standing Height + Speed + Stability = Authentic iLevel (like no other)

iLevel® power seat elevation technology features a revolutionary new power chair design for a better experience. iLevel improves daily functionality and quality of life by allowing operation of the power chair with the seat elevated while at walking speed.

# **FEATURES**

- 30.4 cm (12") of power adjustable seat height while at walking speed; maximum of 25.4 cm (10") height on the Q6 Edge® HD
- iLevel's patented Extra Stability Technology® enhances safety for transfers, reach, and numerous other activities of daily living
- Retrofittable on the Edge® 3 and Q6 Edge® 2.0
- Front and rear LED fender lights come standard with iLevel®





New sleek back design

Optional 30.4 cm (12") iLevel® power seat elevation at 5.6 km/h

Toolless removal of rear shroud with battery access door

> Independent SRS (Smooth Ride Suspension)

NARROWEST POWER BASE IN AUSTRALIA

51.9 CM WIDTH WITH 31.5 CM (12.5") DRIVE WHEELS

SPECS

seat height.

MAXIMUM SPEED: 10 km/h (5.6 km/h when elevated with iLevel®)

The Edge 3 Stretto™ is the narrowest, most manoeuvrable power base in

Australia! Aside from its narrow width, it is equipped with independent

SRS (Smooth Ride Suspension), making it ideal for children, teenagers,

and small adults. The Stretto is available with 31.8 cm (12.5") or 35.6 cm

(14") drive wheels and optional 30.4 cm (12") of iLevel® power adjustable

EDGE 3 STRETTO™


MID-WHEEL DRIVE

WEIGHT CAPACITY: 136 kg

TURNING RADIUS: 53.3 cm

POWER SEAT ELEVATION: Optional 30.4 cm iLevel

# COLOURS

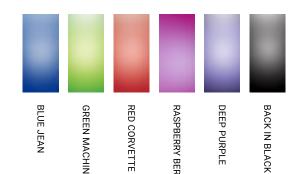

Front and rear LED fender lights

Q-Logic 3 Advanced Drive Controls with

Bluetooth® and USB

technology

Stealth Products® cushion

featuring exclusive Coolcore®

Fun Fact...

Stretto is the Italian word for *narrow*.

Independent SRS (Smooth Ride Suspension)

The Stretto is the ideal paediatric power chair solution

Toolless access to electronics







"iLEVEL HAS GIVEN ME MORE FREEDOM. I LOVE GOING FAST UP HIGH AT 7.24 KM/H."

# QUANTUM EXPERIENCES PROUDLY ASSEMBLED IN AUSTRALIA

As an industry innovator, we are creating a product that offers confidence and stability when elevated at the speed of life. Our Melbourne facility has a robust team of experienced technicians who custom build each Edge 3 specifically to the individual consumer's needs.

# EDGE® 3

# MID-WHEEL DRIVE

The Edge® 3 offers upgraded SRS (Smooth Ride Suspension), which provides a smooth, comfortable ride and advanced stability. Optional iLevel® is available at 7.24 km/h (4.5 mph) with up to 30.4 cm (12") of power adjustable seat height.

# SPECS

MAXIMUM SPEED: 10 km/h (7.24 km/h when elevated with iLevel®)


WEIGHT CAPACITY: 136 kg (See Q6 Edge® HD for increased capacity)

TURNING RADIUS: 52 cm

POWER SEAT ELEVATION: Optional 30.4 cm iLevel

# COLOURS

\* See colour guide for a complete list of options.

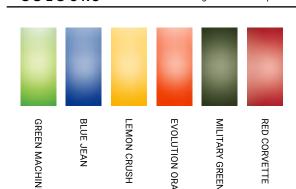

Upgraded SRS (Smooth Ride Suspension)

USB charger (standard)


Featuring Extra Stability
Technology® for safer
transfers, reach, and
other activities of
daily living

QUANTUM®

#1 FOR REHAB POWER

Front and rear LED fender lights for seeing and being seen (standard) Drive wheel colour accents (standard)



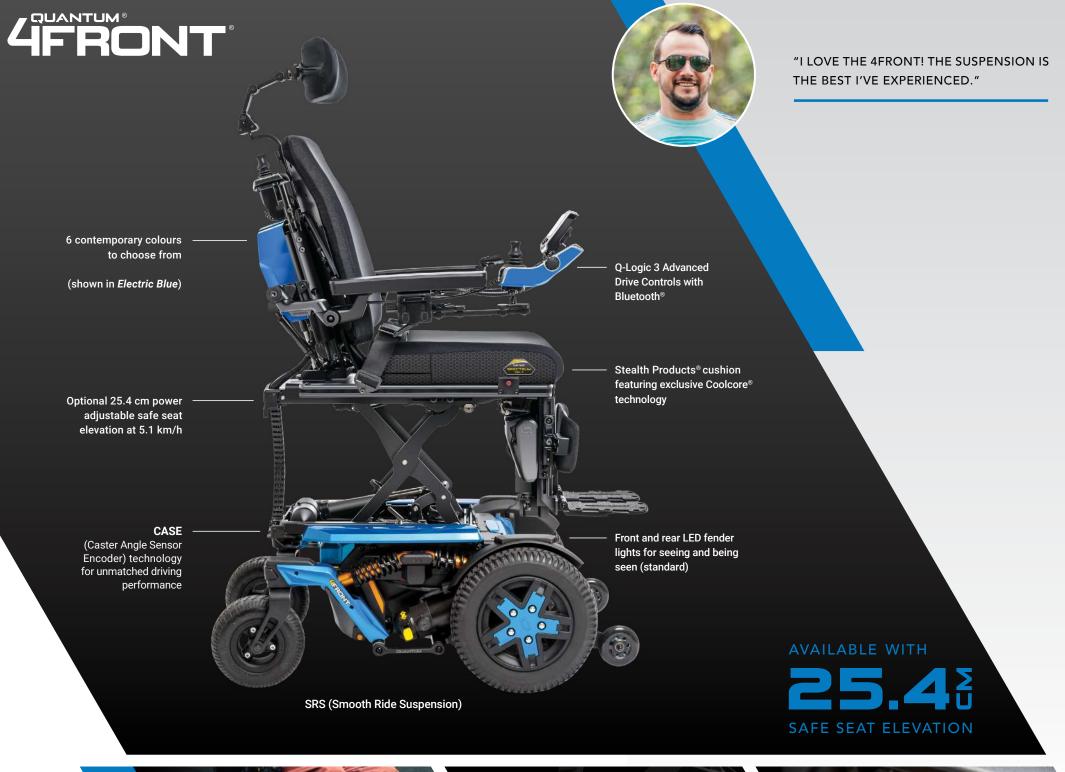

# 4FRONT®

FRONT-WHEEL DRIVE

# Featuring SRS (Smooth Ride Suspension)

The 4Front® is a front-wheel drive power chair that features automotivegrade suspension with unprecedented comfort and ride quality.

# **SPECS**


MAXIMUM SPEED: Up to 10 km/h (5.1 km/h when elevated)

WEIGHT CAPACITY: 136 kg

TURNING RADIUS: 62.25 cm

POWER SEAT ELEVATION: Optional 25.4 cm Safe Seat Elevation

# COLOURS

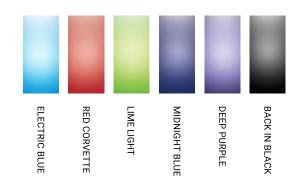









Q6 EDGE® 2.0

MID-WHEEL DRIVE

SPECS

MAXIMUM SPEED: 10 km/h (5.6 km/h when elevated)

WEIGHT CAPACITY: 136 kg

TURNING RADIUS: 52 cm (with 15.2 cm castors)

POWER SEAT ELEVATION: Optional 30.4 cm iLevel®

COLOUR

\* See colour guide for a full list of options.







Q6 EDGE® Z

MID-WHEEL DRIVE

SPECS

MAXIMUM SPEED: 10 km/h (5.6 km/h when elevated)

WEIGHT CAPACITY: 136 kg

TURNING RADIUS: 50 cm

POWER SEAT ELEVATION: Optional 30.4 cm iLevel®

COLOURS



J623 2.0

MID-WHEEL DRIVE

SPECS

MAXIMUM SPEED: 10 km/h

WEIGHT CAPACITY: 136 kg

TURNING RADIUS: 52.1 cm

COLOURS

\* See colour guide for a full list of options.





Q6 EDGE® HD

BARIATRIC MID-WHEEL DRIVE

SPECS

MAXIMUM SPEED: 7.5 km/h (5.6 km/h when elevated)

WEIGHT CAPACITY: 204.1 (181.4 kg with iLevel®)

TURNING RADIUS: 50 cm

POWER SEAT ELEVATION: Optional 25.4 cm iLevel®

COLOURS





R44 RIVAL®

**REAR-WHEEL DRIVE** 

SPECS

MAXIMUM SPEED: 9.6 km/h

WEIGHT CAPACITY: 136 kg

TURNING RADIUS: 90.6 cm

COLOURS

\* See colour guide for a full list of options.



Q1450

**BARIATRIC FRONT-WHEEL DRIVE** 

SPECS

MAXIMUM SPEED: 7.2 km/h

WEIGHT CAPACITY: 272 kg

TURNING RADIUS: 67.3 cm

COLOURS



KOZMO® PORTABLE

PAEDIATRIC REAR-WHEEL DRIVE

SPECS

MAXIMUM SPEED: 5.7 km/h

WEIGHT CAPACITY: 90.7 kg

TURNING RADIUS: 64.6 cm

COLOURS





**QUANTUM®** #1 FOR REHAB POWER

# TRU-BALANCE® 3

# **POWER POSITIONING SYSTEMS**

# **FEATURES**

- Fine adjustments can be made even while client is in the chair
- Seat width can be ordered independently of the back width at no additional charge
- Power tilt, power recline, and iLevel® are retrofittable

# **POWER RECLINE**

85° - 175° of recline	Available in matching back shroud colours
136 kg weight capacity	No increase in seat-to-floor height
Available with power tilt and power adjustable seat lift	Standard shear reduction and cantilever armrests

# **POWER TILT**

50° of tilt	Available in matching back shroud colours		
Lower seat-to-floor height	Maintains center of gravity shift		

Unmatched adjustment and durability

# **HEAVY DUTY POWER TILT & RECLINE**

204 kg weight capacity	Available in matching back shroud colours		
Up to 50° of tilt	85° - 175° of recline		



# Synergy® Seating **Contoured Seating** #1 FOR REHAB POWER

# SEATING

# AND POSITIONING



# SYNERGY® SEATING

Easily adjustable, this versatile design can grow with the user

Several width/depth adjustment packages ranging from 10" to 32"

Accepts a wide variety of after-market seats and backs

Heavy duty option provides a weight capacity of 272 kg and seating up to  $32^{\circ}$  x  $28^{\circ}$  (model specific)



# **COMING SOON**

# CONTOURED SEATING

Black vinyl captain's seat with optional armrests

Compatible with TRU-Balance 3 power articulating foot platform

Available in high back, medium back, solid seat pan, limited recline, and depth-adjustable



# STEALTH PRODUCTS®

**CUSHIONS** 

Stealth Products® Cushions feature patented Coolcore® technology, a chemical free, four-way stretch material that offers moisture wicking, moisture transportation, and moisture evaporation properties. Our cushions also use a colour guide that identifies the suggested function family by a unique patch.

### **FEATURES**

- Chemical free
- Regulates temperature
- Manages and distributes your body heat













QUANTUM®
#1 FOR REHAB POWER

Convenient built-in side/front pouches

Constructed with FRG Sunmate flame-retardant material

Easy-pull rubber loop zipper handles







STEALTH® **GLACIAL SP** 



STEALTH® ZEN SPP

Four-way stretch poly liner

ABS or ACM solid seat insert

Urethane liner



STEALTH® PREMIERE P

Four-way stretch poly liner

ABS or ACM solid seat insert

Urethane liner



STEALTH® SOLUTION® SPP



STEALTH® SPECTRUM® FOAM OR GEL SPP



TRU-COMFORT 2 SPP

Four-way stretch poly liner



STEALTH® **ESSENCE SPP** 

Small contour wedge (0.5")

Four-way stretch poly liner

Medium contour wedge (1.5")

	OPTIONAL:	OPTIONAL:	OPTIONAL:		OPTIONAL:	OPTIONAL:	OPTIONAL:
		Standard visco topper from FRG Sunmate	Standard visco topper from FRG Sunmate			Custom sizes available	
		Made from FRG SunMate flame retardant material	Made from FRG SunMate flame retardant material	Urethane liner for incontinence	Urethane liner for incontinence	Convenient side pouches on both sides of cushion	
Moisture-wicking and Soft gel top layer for skin easy-to-clean Stealth-Tek protection incontinence cover		FRG XXSoft foam topper for skin protection	Built-in positioning platform for off-loading	Twin-cell, aqueous gel insert under at risk bony areas	Multi-layer, viscoelastic foam insert or twin-cell, lightweight gel insert under at risk bony areas	Easy-pull rubber loop zipper handles	Solid seat insert included
	Coolcore® technology cover	Coolcore® technology cover	Coolcore® technology cover	Coolcore® technology cover	Coolcore® technology cover	Coolcore® technology cover	Coolcore® technology cover
Coccyx and sacral relief	Coccyx and sacral relief	Coccyx and sacral relief	Coccyx and sacral relief				Coccyx and sacral relief
Contoured, high-density molded foam base	Contoured, high-density molded foam base	Slightly contoured, high-density molded foam base	Lateral and pre-ischial shelf contoured, high-density molded foam base	Lateral and pre-ischial shelf contoured, high-density molded foam base	Lateral and pre-ischial shelf contoured, high-density molded foam base	Deep contoured, high-density molded foam base	Lateral and pre-ischial shelf contoured, high-density molded foam base
10"W x 10"D - 24"W x 24"D	10"W x 10"D - 24"W x 24"D	10"-24" W x 10"-24" D	16"W x 16"D - 22"W x 22"D	10"W x 10"D - 24"W x 24"D	Foam: 16"W x 15"D - 32"W x 28"D Gel: 14"W x 18"D - 24"W x 26"D	16"W x 16"D - 22"W x 24"D	16"W x 16"D - 22"W x 24"D



# ADI Q-Back 2.0

Four-way stretch poly liner

ABS or ACM solid seat insert

Urethane liner

The ADI Q-Back 2.0 from Stealth Products® is available on TRU-Balance® 3 Seating Systems and is offered in three configurations:

- Active Contour ADI Q-Back 2.0
- Deep Contour ADI Q-Back 2.0
- Extra Deep Contour ADI Q-Back 2.0
- Coolcore® cover comes standard

# Sizes for static/power tilt/power recline seating:

- Support heights: 13", 16", 18", 20", 22"
- Widths: 14", 15", 16", 17", 18", 19", 20", 22"
- Additional back sizes are available with special order

**Note:** 3" of back support height adjustment built into the TRU-Balance® 3 back spine.

# UniLink™ Hardware System

• New system from Stealth Products® for mounting thigh guides and elbow blocks onto the TRU-Balance® 3 Seating System

Four-way stretch poly liner

- Elbow blocks in 3"x10", 3"x3", and 4"x4" pad sizes are made with patented Coolcore® technology cover material
- Thigh guides are available in 4"x4", 4"x8", and 4"x12" pad sizes

# THIGH GUIDE OPTIONS





# TRU-Comfort 2 Back

- · Adjustable positioning and pressure management
- Compatible with both TRU-Balance® 3 Power Positioning Systems and the static seat frame system
- Four-way stretch, recovery, fluid-resistant, breathable cover



# INTERACTIVE ASSIST

# Reduce service visits with our latest innovative app.


- SCREEN MIRRORING
- REAL-TIME MONITORING
- EFFICIENCY IMPROVEMENTS
- TROUBLESHOOTING
- BATTERY MANAGEMENT
- DIAGNOSTIC HISTORY • BASIC PROGRAMMING
- REMOTE ACCESS

- CHANGE LOG HISTORY
  - SAVING XPA FILES
  - MOTOR USAGE HISTORY
  - TEXT CHATTING
  - PROGRAMMING IR MENU
  - ACCESS TO BLUETOOTH FEATURES
  - DEALER-TO-DEALER PROGRAMMING ACCESS

Q-Logic 3 now features Interactive Assist. Interactive Assist from Quantum Rehab® is an app which provides a direct real-time connection from a power chair's electronics to a provider's programming station. The connection provides the technician with full system and diagnostics information along with tools, such as a real-time mirror image of the power chair's electronics display, to eliminate wasteful trips.

**Download Interactive Assist today!** 



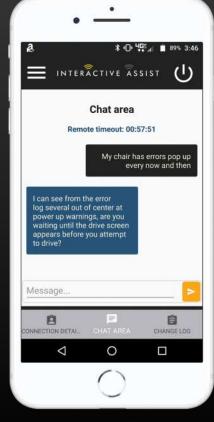



"MY EDGE 3 GIVES ME THE OPPORTUNITY TO USE ADVANCED TECHNOLOGY THAT GIVES ME SO MUCH INDEPENDENCE."



Q-LOGIC 3 EX JOYSTICK CONTROLLER



home screen





ENHANCED DISPLAY

Q-Logic 3 Specialty Controls



Q-LOGIC 3

### ADVANCED DRIVE CONTROLS

The Q-Logic 3 Advanced Drive Controls System from Quantum Rehab® is like no other power wheelchair electronics system. Designed from clinician and consumer feedback, Q-Logic 3 is user-friendly and highly customisable.

### **FEATURES**

- Wireless Bluetooth® programming through Windows devices
- Bluetooth access to computer, tablet, and phone functions
- iAccess programmable push button and toggle options for customised seating access
- Clinic Mode:
  - Quickly set up multiple input devices in client evaluations
  - · Automatically assigned drive modes to input devices

11:25am

Indoor Moderate



Q-Logic 3



drive screen

# **ENHANCED DISPLAY**

The optional Q-Logic 3 EX Enhanced Display serves as a separate screen used with specialty drive controls to display system data. With the Enhanced Display's built-in environmental controls, Q-Logic 3 can operate a wide range of appliances that use infrared (IR) remote controls, such as televisions, and serves as a Bluetooth® computer mouse emulator.









JOYSTICK OPTIONS SIP-N-PUFF MULTI-SWITCH













Stealth® i-Drive® Sip-N-Puff

Stealth® i-Drive® 4 Stealth Pro, 3-Switch Head Array

Stealth® i-Drive® Tri-Array, 3-Switch Head Array

Stealth® i-Drive® Curved Array

TRAY ARRAYS

# **EXPANDABLE (Q-LOGIC 3)**



Q-Logic 3 EX Joystick



Q-Logic 3 EX Enhanced Display



Q-Logic 3 Stand Alone Joystick



Q-Logic 3 iAccess Programmable Switchbox



Q-Logic 3 Specialty Interface Module



SINGLE SWITCH

Stealth® Egg Switch



Mo-Vis™ Twister Basic Button



Mo-Vis<sup>™</sup> Twister Satellite Button



Stealth® i-Drive® Eclipse Tray with Four Capacitive Proximity Sensors



Stealth® i-Drive® Eclipse Tray with Four Fiber Optic Proximity Sensors

# PROPORTIONAL DRIVE CONTROLS



Mo-Vis™ All-Round Joystick



Stealth® Precision Mini Proportional Joystick



Mo-Vis™ All-Round Light Joystick



Stealth® Mushroom Joystick



Mo-Vis™ Multi Joystick



Mo-Vis™ Micro Joystick



Stealth® i-Drive® Control Harness



Front and Rear LED Fender Lights



Personal Item Hook



Hydration System



Phone Holder



Transfer Bars (Pair)



Rearview Mirror



Glove Box

MISC. ACCESSORIES

Quantum® Backpack



# SUPPLEMENTARY

### RESOURCES

At Quantum Rehab®, we are dedicated to providing excellent support services to our customers to ensure you get the most out of your Quantum product. This includes the latest marketing literature and other important educational resources. Contact us or your nearest Authorised Quantum Rehab Provider to learn more.

### **INTRODUCING PQ-365**

PQ-365 includes all the Pride Mobility® and Quantum Rehab® resources needed by providers and clinicians. The free, easy-to-use app helps busy professionals stay up to date with the latest information from Pride and Quantum and access all our content in one quick and convenient place. Features include serial number and order number lookups, a quick-search option, 360-degree product views, fillable order forms, push notifications, consumer testimonials, news, events, specifications sheets, brochures, owner's manuals, installation instructions, basic operating instructions, and more. The app also includes access to the company's Learning Management System (LMS), which provides online and on-demand educational training for professionals.

Available now for Apple and Android devices.











# 26 Q6 EDGE® HD Q6 EDGE® R44 **EDGE®** J623 2.0 STRETTO 4FRONT® EDGE RIVAL® 2.0 Deep Purple Sugar Plum Midnight Blue Blue Jean Electric Blue Military Green Night Train Back in Black **Grey Street** • White Satin • Green Machine • • Lemon Crush Evolution Orange Red Corvette Cherry Bomb Raspberry Beret Pretty in Pink • Lime Light

Shown in Raspberry Beret

# SERVICE

# AND SUPPORT



### STEALTHPRODUCTS.COM

Stealth Products® offers a range of versatile products, including head and body support/positioning, joystick mounts, trays, seating, cushions, and paediatric mobility products.



### LIFEATILEVEL.COM

Learn about the real life adventures and achievements of iLevel® users or share your own. The site features biographies, videos, and a blog.

### LIFE.QUANTUMREHAB.COM

Your one-stop resource for articles and info on Quantum's® life-changing power wheelchairs and complex rehab technologies.










# SOCIAL MEDIA

Stay up to date on the latest Quantum® events and product releases through Twitter, Facebook, YouTube, and Instagram.

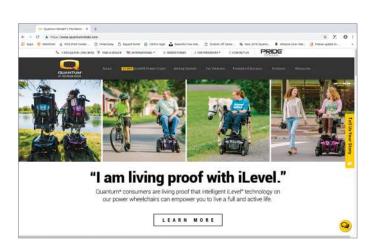

### QUANTUMREHAB.COM

Quantum's official website is the ultimate source for everything you need to know about Quantum® products. Highlights include detailed product specifications, order forms, videos, user testimonials, and much more.

### QUANTUM REHAB PROFESSIONAL APP

The Quantum Rehab® Professional App is available for Apple and Android smart phones and tablets. Look for the app in the Apple Store or Google Play using your smart device to download the app directly.

# SALES & TECHNICAL SUPPORT

For more information, contact your local Authorised Quantum Provider.

Stamp/Quantum® Provider Details

<sup>\*</sup> Additional J623 2.0 colours include: Carbon Fibre, Emerald Green, Orange Crush, Sapphire Blue, Silver Bullet, & Yellow Blaze. Contact your nearest Authorised Quantum Provider to learn more.



AVAILABLE WITH 30.4 CM iLEVEL®

AVAILABLE WITH 30.4 CM iLEVEL®

AVAILABLE WITH 30.4 CM iLEVEL®

AVAILABLE WITH 25.4 CM ELEVATION





#1 FUR REHAB PUWER	EDGE 3 STRETTO®	EDGE® 3	Q6 EDGE® 2.0	4FRONT®	
DRIVE CONFIGURATION	Mid-wheel Drive	Mid-wheel Drive	Mid-wheel Drive	Front-wheel Drive	
DRIVE WHEELS	35.5 cm (14") or 31.7 cm (12.5") Pneumatic or Flat Free Tyres	35.5 cm (14") Flat Free/ 35.5 cm (14") Pneumatic	35.5 cm (14") Flat Free/ 35.5 cm (14") Pneumatic	35.5 cm (14") Flat Free/ 35.5 cm (14") Pneumatic	
CASTOR WHEELS	Front: 12.7 cm (5") Flat Free	Front: 15.2 cm (6") Solid	<b>Front:</b> 15.2 cm (6") Solid <sup>1</sup> 20.3 cm (8") Solid <sup>2</sup>	- <b>Rear:</b> 22.8 cm (9") Flat Free	
OAGTOR IIIIEEEG	Rear: 15.2 cm (6") Flat Free	<b>Rear:</b> 15.2 cm (6") Solid	<b>Rear:</b> 15.2 cm (6") Solid <sup>1</sup> 20.3 cm (8") Solid <sup>2</sup>		
ANTI-TIP WHEELS	N/A	N/A	N/A	Front: 10.1 cm (4")	
MAXIMUM SPEED <sup>4</sup>	Up to 10 km/h (6.2 mph) (5.6 km/h (3.5 mph) w/iLevel) <sup>13</sup>	Up to 10 km/h (6.2 mph)	Up to 10 km/h (6.2 mph)	Up to 10 km/h (6.2 mph)	
GROUND CLEARANCE	7.3 cm (2.9")	6.8 cm (2.7") (battery tray)	6.8 cm (2.7") (battery tray)	7.6 cm (3") (motor pivot)	
TURNING RADIUS <sup>5</sup>	53.3 cm (21")	52 cm (20.5")	52 cm (20.5") w/15.2 cm (6") castors 59 cm (23.2") w/20.3 cm (8") castors	62.2 cm (24.5")	
OVERALL LENGTH <sup>11</sup>	82.5 cm (32.5")	81.8 cm (32.2") (without foot riggings)	90.1 cm (35.5") w/15.2 cm (6") castors 102.6 cm (40.4") w/20.3 cm (8") castors	116 cm (45.7")	
55.2 cm (21.75") (35.5 cm (14") drive wheels)  51.9 cm (20.47") (31.5 cm (12.5") drive wheels)		61.9 cm (24.4")	61.9 cm (24.4")	60.9 cm (24")	
DASE WIDTH	51.9 cm (20.47") (31.5 cm (12.5") drive wheels)	01.5011(21.17)	01.3 011 (24.47)	55.5 SIN (E 1 )	
SEAT SIZES	, , ,	TRU-Balance® 3	TRU-Balance® 3	TRU-Balance® 3	
	(12.5") drive wheels)	, ,		, ,	
SEAT SIZES	(12.5") drive wheels)  TRU-Balance® 3	TRU-Balance® 3	TRU-Balance® 3	TRU-Balance® 3	
SEAT SIZES SYNERGY®	(12.5") drive wheels)  TRU-Balance® 3  N/A	TRU-Balance® 3	TRU-Balance® 3 W: 12-22" D: 12-22"	TRU-Balance® 3	
SEAT SIZES SYNERGY® STATIC	(12.5") drive wheels)  TRU-Balance® 3  N/A  W: 12-20" D: 12-20"	TRU-Balance® 3 N/A W: 12-22" D: 10-22"	TRU-Balance® 3 W: 12-22" D: 12-22" W: 12-22" D: 12-22"	TRU-Balance® 3 N/A W: 12-22" D: 12-22"	
SEAT SIZES SYNERGY® STATIC POWER TILT	(12.5") drive wheels)  TRU-Balance® 3  N/A  W: 12-20" D: 12-20"  W: 12-20" D: 12-20"  W: 12-20" D: 12-20"	TRU-Balance® 3  N/A  W: 12-22" <b>D</b> : 10-22"  W: 12-22" <b>D</b> : 12-22"  W: 12-22" <b>D</b> : 12-22"	TRU-Balance® 3  W: 12-22" D: 12-22"  W: 12-22" D: 12-22"  W: 12-22" D: 12-22"  W: 12-22" D: 12-22"	TRU-Balance® 3  N/A  W: 12-22" D: 12-22"  W: 12-22" D: 12-22"  W: 12-22" D: 12-22"  TRU-Balance® 3	
SEAT SIZES SYNERGY® STATIC POWER TILT LIFT & TILT SEAT-TO-FLOOR HEIGHTS 16,17 SYNERGY®/STATIC	(12.5") drive wheels)  TRU-Balance® 3  N/A  W: 12-20" D: 12-20"  W: 12-20" D: 12-20"	TRU-Balance® 3  N/A  W: 12-22" <b>D</b> : 10-22"  W: 12-22" <b>D</b> : 12-22"  W: 12-22" <b>D</b> : 12-22"  TRU-Balance® 3  16.5-18.5"	TRU-Balance® 3  W: 12-22" D: 12-22"  W: 12-22" D: 12-22"  W: 12-22" D: 12-22"  TRU-Balance® 3  16.5-18.5"	TRU-Balance® 3  N/A  W: 12-22" D: 12-22"  W: 12-22" D: 12-22"  W: 12-22" D: 12-22"  TRU-Balance® 3  15.2-17.2"3	
SEAT SIZES SYNERGY® STATIC POWER TILT LIFT & TILT SEAT-TO-FLOOR HEIGHTS 16.17 SYNERGY®/STATIC POWER TILT®	(12.5") drive wheels)  TRU-Balance® 3  N/A  W: 12-20" <b>D</b> : 12-20"  W: 12-20" <b>D</b> : 12-20"  W: 12-20" <b>D</b> : 12-20"  TRU-Balance® 3  16"-18" (Static only)  17.5-19.5"	TRU-Balance® 3  N/A  W: 12-22" <b>D</b> : 10-22"  W: 12-22" <b>D</b> : 12-22"  W: 12-22" <b>D</b> : 12-22"  TRU-Balance® 3  16.5-18.5"  17.5-19.5"	TRU-Balance® 3  W: 12-22" D: 12-22"  W: 12-22" D: 12-22"  W: 12-22" D: 12-22"  W: 12-22" D: 12-22"  TRU-Balance® 3  16.5-18.5"  17.5-19.5"	TRU-Balance® 3  N/A  W: 12-22" D: 12-22"  W: 12-22" D: 12-22"  W: 12-22" D: 12-22"  TRU-Balance® 3  15.2-17.2" <sup>3</sup> 17.5-19.5"	
SEAT SIZES SYNERGY® STATIC POWER TILT LIFT & TILT SEAT-TO-FLOOR HEIGHTS 16.17 SYNERGY®/STATIC POWER TILT® LIFT & TILT®	(12.5") drive wheels)  TRU-Balance® 3  N/A  W: 12-20" D: 12-20"  W: 12-20" D: 12-20"  W: 12-20" D: 12-20"  TRU-Balance® 3  16"-18" (Static only)  17.5-19.5"	TRU-Balance® 3  N/A  W: 12-22" D: 10-22"  W: 12-22" D: 12-22"  W: 12-22" D: 12-22"  TRU-Balance® 3  16.5-18.5"  17.5-19.5"	TRU-Balance® 3  W: 12-22" D: 12-22"  W: 12-22" D: 12-22"  W: 12-22" D: 12-22"  TRU-Balance® 3  16.5-18.5"  17.5-19.5"	TRU-Balance® 3  N/A  W: 12-22" D: 12-22"  W: 12-22" D: 12-22"  W: 12-22" D: 12-22"  TRU-Balance® 3  15.2-17.2"3  17.5-19.5"	
SEAT SIZES SYNERGY® STATIC POWER TILT LIFT & TILT SEAT-TO-FLOOR HEIGHTS 16,17 SYNERGY®/STATIC POWER TILT® LIFT & TILT® LIFT & TILT® MANUAL RECLINE	(12.5") drive wheels)  TRU-Balance® 3  N/A  W: 12-20" D: 12-20"  W: 12-20" D: 12-20"  TRU-Balance® 3  16"-18" (Static only)  17.5-19.5"  Yes	TRU-Balance® 3  N/A  W: 12-22" <b>D</b> : 10-22"  W: 12-22" <b>D</b> : 12-22"  W: 12-22" <b>D</b> : 12-22"  TRU-Balance® 3  16.5-18.5"  17.5-19.5"  Yes	TRU-Balance® 3  W: 12-22" D: 12-22"  W: 12-22" D: 12-22"  W: 12-22" D: 12-22"  W: 12-22" D: 12-22"  TRU-Balance® 3  16.5-18.5"  17.5-19.5"  Yes	TRU-Balance® 3  N/A  W: 12-22" D: 12-22"  W: 12-22" D: 12-22"  W: 12-22" D: 12-22"  TRU-Balance® 3  15.2-17.2"3  17.5-19.5"  Yes	
SEAT SIZES SYNERGY® STATIC POWER TILT LIFT & TILT SEAT-TO-FLOOR HEIGHTS 16.17 SYNERGY®/STATIC POWER TILT® LIFT & TILT®	(12.5") drive wheels)  TRU-Balance® 3  N/A  W: 12-20" D: 12-20"  W: 12-20" D: 12-20"  TRU-Balance® 3  16"-18" (Static only)  17.5-19.5"  Yes  40AH (2) or NF-22 (55AH) (2)	TRU-Balance® 3  N/A  W: 12-22" D: 10-22"  W: 12-22" D: 12-22"  W: 12-22" D: 12-22"  TRU-Balance® 3  16.5-18.5"  17.5-19.5"	TRU-Balance® 3  W: 12-22" D: 12-22"  W: 12-22" D: 12-22"  W: 12-22" D: 12-22"  TRU-Balance® 3  16.5-18.5"  17.5-19.5"	TRU-Balance® 3  N/A  W: 12-22" D: 12-22"  W: 12-22" D: 12-22"  W: 12-22" D: 12-22"  TRU-Balance® 3  15.2-17.2"3  17.5-19.5"	
SEAT SIZES SYNERGY® STATIC POWER TILT LIFT & TILT SEAT-TO-FLOOR HEIGHTS 16,17 SYNERGY®/STATIC POWER TILT® LIFT & TILT® LIFT & TILT® MANUAL RECLINE	(12.5") drive wheels)  TRU-Balance® 3  N/A  W: 12-20" D: 12-20"  W: 12-20" D: 12-20"  TRU-Balance® 3  16"-18" (Static only)  17.5-19.5"  Yes	TRU-Balance® 3  N/A  W: 12-22" <b>D</b> : 10-22"  W: 12-22" <b>D</b> : 12-22"  W: 12-22" <b>D</b> : 12-22"  TRU-Balance® 3  16.5-18.5"  17.5-19.5"  Yes	TRU-Balance® 3  W: 12-22" D: 12-22"  W: 12-22" D: 12-22"  W: 12-22" D: 12-22"  W: 12-22" D: 12-22"  TRU-Balance® 3  16.5-18.5"  17.5-19.5"  Yes	TRU-Balance® 3  N/A  W: 12-22" D: 12-22"  W: 12-22" D: 12-22"  W: 12-22" D: 12-22"  TRU-Balance® 3  15.2-17.2"3  17.5-19.5"  Yes	
SEAT SIZES SYNERGY® STATIC POWER TILT LIFT & TILT SEAT-TO-FLOOR HEIGHTS 16.17 SYNERGY®/STATIC POWER TILT® LIFT & TILT® MANUAL RECLINE BATTERY SIZE®	(12.5") drive wheels)  TRU-Balance® 3  N/A  W: 12-20" D: 12-20"  W: 12-20" D: 12-20"  TRU-Balance® 3  16"-18" (Static only)  17.5-19.5"  Yes  40AH (2) or NF-22 (55AH) (2)  40 AH: 14 kgs (31 lbs. each)	TRU-Balance® 3  N/A  W: 12-22" <b>D</b> : 10-22"  W: 12-22" <b>D</b> : 12-22"  W: 12-22" <b>D</b> : 12-22"  TRU-Balance® 3  16.5-18.5"  17.5-19.5"  Yes  NF-22 (55AH) (2)	TRU-Balance® 3  W: 12-22" D: 12-22"  W: 12-22" D: 12-22"  W: 12-22" D: 12-22"  W: 12-22" D: 12-22"  TRU-Balance® 3  16.5-18.5"  17.5-19.5"  Yes  NF-22 (55AH) (2)	TRU-Balance® 3  N/A  W: 12-22" <b>D</b> : 12-22"  W: 12-22" <b>D</b> : 12-22"  TRU-Balance® 3  15.2-17.2"³  17.5-19.5"  Yes  Group 24: (74AH) (2)	
SEAT SIZES SYNERGY® STATIC POWER TILT LIFT & TILT SEAT-TO-FLOOR HEIGHTS 16,17 SYNERGY®/STATIC POWER TILT® LIFT & TILT® MANUAL RECLINE BATTERY SIZE® BATTERY WEIGHT14	(12.5") drive wheels)  TRU-Balance® 3  N/A  W: 12-20" D: 12-20"  W: 12-20" D: 12-20"  W: 12-20" D: 12-20"  TRU-Balance® 3  16"-18" (Static only)  17.5-19.5"  Yes  40AH (2) or NF-22 (55AH) (2)  40 AH: 14 kgs (31 lbs. each)  55 AH: 18.1 kgs (40 lbs. each)	TRU-Balance® 3  N/A  W: 12-22" <b>D</b> : 10-22"  W: 12-22" <b>D</b> : 12-22"  W: 12-22" <b>D</b> : 12-22"  TRU-Balance® 3  16.5-18.5"  17.5-19.5"  Yes  NF-22 (55AH) (2)  17.2 kgs (37.9 lbs.) (each)	TRU-Balance® 3  W: 12-22" D: 12-22"  W: 12-22" D: 12-22"  W: 12-22" D: 12-22"  TRU-Balance® 3  16.5-18.5"  17.5-19.5"  Yes  NF-22 (55AH) (2)  17.2 kgs (37.9 lbs.) (each)	TRU-Balance® 3  N/A  W: 12-22" D: 12-22"  W: 12-22" D: 12-22"  W: 12-22" D: 12-22"  TRU-Balance® 3  15.2-17.2"3  17.5-19.5"  Yes  Group 24: (74AH) (2)  23.5 kgs (51.8 lbs.) (each)	
SEAT SIZES SYNERGY® STATIC POWER TILT LIFT & TILT SEAT-TO-FLOOR HEIGHTS 16,17 SYNERGY®/STATIC POWER TILT® LIFT & TILT® MANUAL RECLINE BATTERY SIZE® BATTERY WEIGHT 14  AVAILABLE ELECTRONICS 12	(12.5") drive wheels)  TRU-Balance® 3  N/A  W: 12-20" D: 12-20"  W: 12-20" D: 12-20"  W: 12-20" D: 12-20"  TRU-Balance® 3  16"-18" (Static only)  17.5-19.5"  Yes  40AH (2) or NF-22 (55AH) (2)  40 AH: 14 kgs (31 lbs. each)  55 AH: 18.1 kgs (40 lbs. each)  120A Q-Logic 3 EX	TRU-Balance® 3  N/A  W: 12-22" D: 10-22"  W: 12-22" D: 12-22"  W: 12-22" D: 12-22"  TRU-Balance® 3  16.5-18.5"  17.5-19.5"  Yes  NF-22 (55AH) (2)  17.2 kgs (37.9 lbs.) (each)  120A Q-Logic 3 EX	TRU-Balance® 3  W: 12-22" D: 12-22"  W: 12-22" D: 12-22"  W: 12-22" D: 12-22"  W: 12-22" D: 12-22"  TRU-Balance® 3  16.5-18.5"  17.5-19.5"  Yes  NF-22 (55AH) (2)  17.2 kgs (37.9 lbs.) (each)  120A Q-Logic 3 EX	TRU-Balance® 3  N/A  W: 12-22" <b>D</b> : 12-22"  W: 12-22" <b>D</b> : 12-22"  TRU-Balance® 3  15.2-17.2"3  17.5-19.5"  Yes  Group 24: (74AH) (2)  23.5 kgs (51.8 lbs.) (each)  120A Q-Logic 3 EX	
SEAT SIZES SYNERGY® STATIC POWER TILT LIFT & TILT SEAT-TO-FLOOR HEIGHTS 16,17 SYNERGY®/STATIC POWER TILT® LIFT & TILT® LIFT & TILT® ANAIUAL RECLINE BATTERY SIZE® BATTERY WEIGHT14  AVAILABLE ELECTRONICS12 BATTERY CHARGER	(12.5") drive wheels)  TRU-Balance® 3  N/A  W: 12-20" D: 12-20"  W: 12-20" D: 12-20"  W: 12-20" D: 12-20"  TRU-Balance® 3  16"-18" (Static only)  17.5-19.5"  Yes  40AH (2) or NF-22 (55AH) (2)  40 AH: 14 kgs (31 lbs. each)  55 AH: 18.1 kgs (40 lbs. each)  120A Q-Logic 3 EX  8A Off-board	TRU-Balance® 3  N/A  W: 12-22" D: 10-22"  W: 12-22" D: 12-22"  W: 12-22" D: 12-22"  TRU-Balance® 3  16.5-18.5"  17.5-19.5"  Yes  NF-22 (55AH) (2)  17.2 kgs (37.9 lbs.) (each)  120A Q-Logic 3 EX  8A Off-board	TRU-Balance® 3  W: 12-22" D: 12-22"  W: 12-22" D: 12-22"  W: 12-22" D: 12-22"  TRU-Balance® 3  16.5-18.5"  17.5-19.5"  Yes  NF-22 (55AH) (2)  17.2 kgs (37.9 lbs.) (each)  120A Q-Logic 3 EX  8A Off-board	TRU-Balance® 3  N/A  W: 12-22" D: 12-22"  W: 12-22" D: 12-22"  W: 12-22" D: 12-22"  TRU-Balance® 3  15.2-17.2"3  17.5-19.5"  Yes  Group 24: (74AH) (2)  23.5 kgs (51.8 lbs.) (each)  120A Q-Logic 3 EX  8A Off-board	

<sup>1</sup> Standard. 2 Optional. 3 Not available with Synergy®. 4 Speed varies with user weight, terrain type, battery charge, battery condition and tyre pressure. 5 Without front riggings. 6 Battery capacity may change when adding power positioning systems. 7 Weight capacity is dictated by the selection of seating. 8 All measurements +/- 1/8" with 18" x 18" seat, fixed towers, 136 kg weight capacity and solid tyres (add 3/8" with pneumatic tyres). Depths may increase seat-to-floor height. Contact Customer Service for details. Listed measurements are in their lowest position; does not include elevating system unless specified.  $9\,113\,kg$  weight capacity with TRU-Balance® Systems.  $10\,A$  minimum seating size of  $20"\times20"$  is required with iLevel® on











AVAILABLE WITH

25.4 CM ELEVATION



AVAILABLE WITH

25.4 CM ELEVATION

Q6 EDGE® Z	Q6 EDGE® HD		QUANTUM® 1450	J623 2.0	QUANTUM® R44 RIVAL®	
Mid-wheel Drive	Mid-wheel Drive		Front-wheel Drive	Mid-wheel Drive	Rear-wheel Drive	
35.5 cm (14") Flat Free/ 35.5 cm (14") Pneumatic	35.5 cm (14") Flat Free/ 35.5 cm (14") Pneumatic		35.5 cm (14") Flat Free	35.5 cm (14") Flat Free/ 35.5 cm (14") Pneumatic	35.5 cm (14") Flat Free/ 35.5 cm (14") Pneumatic	
<b>Front:</b> 15.2 cm (6") Solid	Front: 15.2 cm (6") Solid Rear: 15.2 cm (6") Solid		22.8 cm (9") Solid	Front: 15.2 cm (6") Solid <sup>1</sup> 20.3 cm (8") Solid <sup>2</sup>	Front: 22.8 cm (9") Solid/Pneumatic 25.4 (10") Solid/Pneumatic	
<b>Rear:</b> 15.2 cm (6") Solid			22.0 0111 (3 ) 30110	<b>Rear:</b> 15.2 cm (6") Solid <sup>1</sup> 20.3 cm (8") Solid <sup>2</sup>		
N/A	N/A		Front: 10.1 cm (4") Solid	N/A	<b>Rear:</b> 7.6 cm (3") Solid	
Up to 10 km/h (6.2 mph)	Up to 7.5 km/h (4.6 mph)		Up to 7.2 km/h (4.4 mph)	Up to 10 km/h (6.2 mph)	Up to 9.6 km/h (6 mph)	
7.6 cm (3") (frame base)	7.6 cm (3") (frame cen	tre)	5 cm (2") (anti-tips)	6.9 cm (2.7") (battery tray)	8.1 cm (3.2") (frame b	ase)
50 cm (19.7")	50 cm (19.7")		67.3 cm (26.5")	52.1 cm (20.5")	90.6 cm (35.7)	
90.4 cm (35.6")	90.4 cm (35.6")		90.4 cm (35.6")	90.2 cm (35.5") without foot riggings	91.4 cm (36")	
66.5 cm (26.2")	66.5 cm (26.2")		74.1 cm (29.2")	61.9 cm (24.4")	63.5 cm (25")	
TRU-Balance® 3	TB2 HD	TB3 HD	TRU-Balance® 2 HD Tilt	TRU-Balance® 3	TB2 HD	TB3 HD
<b>W</b> : 12-24" <b>D</b> : 10-24"	<b>W</b> : 16-28" <b>D</b> : 20-28"	<b>W</b> : 16-32" <b>D</b> : 20-28"	<b>W</b> : 20-32" <b>D</b> : 20-28"	W: 12-22" D: 12-22"	<b>W</b> : 12-22" <b>D</b> : 10-20"	W: 12-22" D: 12-22"
<b>W</b> : 12-24" <b>D</b> : 12-24"	<b>W</b> : 16-28" <b>D</b> : 20-28"	<b>W</b> : 16-32" <b>D</b> : 20-28"	<b>W</b> : 20-32" <b>D</b> : 20-28"	W: 12-22" D: 12-22"	<b>W:</b> 12-22" <b>D:</b> 10-20"	<b>W:</b> 12-22" <b>D:</b> 12-22"
<b>W</b> : 12-24" <b>D</b> : 14-24"	<b>W</b> : 16-28" <b>D</b> : 20-28"	<b>W</b> : 16-28" <b>D</b> : 20-25"	<b>W</b> : 20-32" <b>D</b> : 20-28"	W: 12-22" D: 12-22"	<b>W</b> : 14-22" <b>D</b> : 14-20"	<b>W</b> : 12-22" <b>D</b> : 12-22"
<b>W</b> : 12-24" <b>D</b> : 12-24"	N/A	<b>W</b> : 16-24" <b>D</b> : 20-25"	N/A	N/A	N/A	<b>W</b> : 12-22" <b>D</b> : 12-22"
TRU-Balance® 3	TB2 HD	TB3 HD	TRU-Balance® 2 HD Tilt	TRU-Balance® 3	TB2 HD	TB3 HD
16.5-18.5"	16.5-18.5"	17.5-19.5"	16.5-18.5"	16.5-18.5"	16.5-18"	16.5-18.1"
17.5-19.5"	18.2-20.9"	18.5-20.5"	18.1-20.1"	17.5-19.5"	16.5-18"	18-19.6"
18-20"	N/A	17.7-19.7"	N/A	17.5-19.5"	N/A	18.5-20.1"
Yes	No		No	Yes	Yes	
Group 24: (74AH) (2)	Group 24: (74AH) (2)		Group 24: (74AH) (2)	NF-22 (55AH) (2)	Group 24: (74AH) (2)	
23.5 kgs (51.8 lbs.) (each)	23.5 kgs (51.8 lbs.) (each)		23.5 kgs (51.8 lbs.) (each)	17.2 kgs (37.9 lbs.) (each)	23.5 kgs (51.8 lbs.) (each)	
120A Q-Logic 3 EX	120A Q-Logic 3 EX		90A VR2 6-key 100A Q-Logic 2 EX	90A VR2 6-key	100A Q-Logic 2 EX	
8A Off-board	8A Off-board		8A Off-board	8A Off-board	8A Off-board	
4-Pole, Mid-Wheel 6®	4-Pole, Mid-Wheel 6®1		4-Pole, Front-wheel Drive	2-Pole Mid-Wheel 6®	4-Pole, Rear-wheel Drive <sup>1</sup>	
136 kgs (300 lbs.)	204.1 kgs (450 lbs.) <sup>15</sup>		272 kgs (600 lbs.)	136 kgs (300 lbs.)	136 kgs (300 lbs.)	
65 kgs (143.3 lbs.)	65 kgs (143.3 lbs.)		88.4 kgs (194.8 lbs.)	62.8 kgs (138.6 lbs.)	58.5 kgs (128.9 lbs.)	

the Q6 Edge® HD; maximum seating size is 24" x 25". 11 Without foot riggings. 12 All Quantum units can be ordered with specialty controls. 13 10 km/h speed available only with Q-Logic 3. 14 Battery weight may vary +/- 1.3 kgs based on manufacturer. 15 Weight capacity is reduced to 181 kgs with iLevel® power adjustable seat height. 16 Seat-to-floor heights are measured with the unit unloaded. Due to the 4Front's suspension system, seat-to-floor heights are approximately 0.6 cm lower when the unit is loaded. 17 41.9 cm seat-to-floor height is not available with power tilt or iLevel®. 49.5 cm seat-to-floor height is not available with static seat or power recline. NOTE: The information contained herein is correct at the time of publication.



#1 FOR REHAB POWER