

4FRONT™

The 4Front™ is a quiet, more responsive front-wheel drive power chair that features automotive-grade suspension with unprecedented comfort and ride quality. The 4Front's CASE (Caster Angle Sensor Encoder) technology and SRS (Smooth Ride Suspension) provide stability for smooth obstacle transitions and unmatched driving performance. The 4Front™ features standard speeds of up to 9.6 km/h (6 mph) in the lowered position, and speeds up to 5.1 km/h (3.2 mph) while elevated at 25.4 cm (10") with our optional safe seat elevation.

QUANTUM® 4FRONT™ shown in Electric Blue with optional *Safe Seat Elevation* and Q-Logic 3 Advanced Drive Controls with Bluetooth®



QUANTUM®

Dedicated to revolutionising the power chair experience.

(TEL) 03 8770 9600

(FAX) 03 9703 2720

www.quantumrehab.com



Specifications

DRIVE WHEELS	35.56 cm (14") flat free
CASTOR WHEELS	Rear: 22.86 cm (9") flat free
ANTI-TIP WHEELS	Front: 10.16 cm (4")
MAXIMUM SPEED ¹	Up to 9.6 km/h (6 mph) Up to 5.1 km/h (3.2 mph) with 25.4 cm (10") lift
GROUND CLEARANCE	6.7 cm (2.6") (motor pivot)
TURNING RADIUS ²	62.25 cm (24.5")
OVERALL LENGTH ²	101.4 cm (39.9")
BASE WIDTH	61.6 cm (24.2")
SEATING SIZES	
Synergy®/TRU-Balance® 3 Static Seat	Width: 30.5-56 cm (12-22") Depth: 30.5-56 cm (12-22")
TRU-Balance® 3 Power Tilt	Width: 30.5-56 cm (12-22") Depth: 30.5-56 cm (12-22")
TRU-Balance® 3 Lift & Tilt	Width: 30.5-56 cm (12-22") Depth: 30.5-56 cm (12-22")
SEAT-TO-FLOOR HEIGHTS ⁹	
Synergy®/TRU-Balance® 3 Static Seat ⁵	38.7-43.8 cm (15.25-17.25")
TRU-Balance® 3 Power Tilt ⁵	44.5-49.5 cm (17.5-19.5")
TRU-Balance® 3 Lift & Tilt ⁵	45-50 cm (17.75-19.75")
MANUAL RECLINE	Yes
BATTERY SIZE ³	Group 24 or Group 34
BATTERY WEIGHT ⁶	23.4 kg (51.5 lbs.) each 17.5 kg (38.5 lbs.) each
AVAILABLE ELECTRONICS ⁷	120A Q-Logic 3 EX
BATTERY CHARGER	8A Off-board
MOTOR PACKAGES	Two-pole, front-wheel drive
WEIGHT CAPACITY ⁴	136 kg (300 lbs.)
BASE WEIGHT	63.5 kg (140 lbs.)

1) Speed varies with user weight, terrain type, battery charge, battery condition and tyre pressure.
2) Without front riggings
3) Battery capacity may change when adding power positioning systems.
4) Weight capacity is dictated by the selection of seating.
5) All measurements +/- 3 mm (1/8") with 46 x 46 cm (18" x 18") seat, fixed towers, 136 kg (300 lbs.) weight capacity and solid tyres (add 10 mm (3/8") with pneumatic tires). Widths may increase seat-to-floor height. Contact Customer Service for details. Listed measurements are in the lowest position; does not include elevating system unless specified.
6) Battery weight may vary +/- 1.4 kg (3 lbs.) based on manufacturer.
7) All Quantum Rehab® units can be ordered with specialty controls.
8) Call for availability
9) Seat-to-floor heights are measured with the unit unloaded. Due to the 4Front's suspension system, seat-to-floor heights are approximately 6 mm (0.25") lower when the unit is loaded.

Standard Features

- CASE (Caster Angle Sensor Encoder) technology
- Unmatched driving performance and ride quality
- Advanced stability for smooth obstacle transitions
- Superior, automotive-grade SRS (Smooth Ride Suspension)
- Q-Logic 3 EX Advanced Drive Controls
- Quiet ride technology




Drive Wheel Colour Accents

Options

- 25.4 cm (10") of power adjustable safe seat elevation at speeds of up to 5.1 km/h (3.2 mph)
- TRU-Comfort 2 Seating
- Trilogy vent tray ⁸
- Compact vent tray ⁸
- TRU-Balance® 3 Seating
- LED lighting package

Colours

SEAT BACKS for units without TRU-Balance® 3 are only available in *Back in Black*.

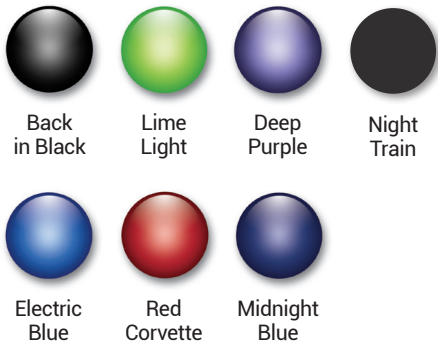

QUANTUM®



20-24 Apollo Drive Hallam, Victoria 3803

ACN 088 609 661 ABN 66 088 609 661

(Tel) 03 8770 9600 • (Fax) 03 9703 2720

www.pridemobility.com.au • www.quantumrehab.com



/quantumrehab



@quantum_rehab

Pride®
 Mobility Products
 Australia Pty. Ltd.
 — LIVE YOUR BEST —

